



EC REP

Shanghai International Holding Corp. GmbH (Europe)
Add: Eiffestrasse 80, 20537 Hamburg, Germany
Tel: +49-40-2513175 Fax: +49-40-255726

Manufacturer: Guangdong JINME Medical Technology Co., Ltd.
Tel: +86-757-88777991
www.appledentalmed.com
2022/12 Version



THE DENTAL
CHAIR & HANDPIECE
COMPANY



AP026

Innovation in totally harmony
European style



Ergonomics concept design

AP026 was designed based on humanized concept, it with 9 position memories, can suit for different users. It created with collision-avoidance system, when chair back hit something it can automatically stop descending, can protect chair or people in time. Its spittoon can be rotated 180°, side box can be rotated 60°, expanding more space for patient or assistant.



The unit box can be rotated 60°



To expand the space for assistant position.

Ergonomic patient positioning



The designed angles of cushion were tested by more than 100 people to let combined movements of the backrest and seat for an ergonomic extension of the spine.

Efficient&Convenient

It is like European style ,which main instrument tray and assistant tray have used touch system, for clear and convenient operate. It have choice for 3 dentists, and each dentist has 3 position memories ,total 9 memories.with “LP” and “R”-Reset,“LP”can make the position to the lowest, and lamp can be automatically off for patient spitting, After that, then push “LP” , it recovers to the position which stayed just now and lamp be on automatically. “R” -Reset is an emergency key, to low down the position to lowest for dealing emergency situation.



Smart-touch panel



AP026 adopts the modern technology of touch screen,to operate the function easily and fastly, 9 momories,with “LP” and “R” functions.

LED reflector sensor lamp



With sensitive induction of this LED lamp,its soft light can makes the patients feel comfortable for eyes during operation time.

Easily adjustable headrest



Unlike the normally screw type to adjust this headrest, this model only need a little strength to push the pole then can get any position you want lying.

Multifunctional foot pedal



Multi-function foot pedal with big size convenient to operate.

The difference is in the details

AP026 have the details which shows more durable -it adopts the USA tubes and Italy valves and Taiwan motor. Also the special Collision-avoidance function can help the dentists or products be protected in time.



European style big side box



With an elegant APPLIEDENTAL logo and Inside with American tubes and Italy valves, can make the chair quality more superior.

Big size backrest



AP026 Adopts European style size of cushion,can suit for different weight of patients,with its strong base can bear max weight 280kg.

Collision-avoidance function



when chair back hit something or meet a certain amount of resistance,it can automatically stop descending,can protect chair or people/products in time.

Ceramic cuspidor



This cuspidor can be rotated for 180°,convenient for the patient to use, and also it is detachable for doing deeply clean work.

UP-Mounted for choice

A nice European design instrument tray has adjustable pull rod and silicone protection pad which can better protect syringe, handpiece and scaler. Screen touch controller provides extremely high sensitivity and ultrafast response which can make your work in cost effective and timely manner.



Tray can be rotated 360°.




Easy and clear touch panel, with three memories, with time and operation data.



up-mounted tool tray tray, the tubes with stepwise stretching and automatically reset function, no dragging feeling.

INSTALLATION DRAWING

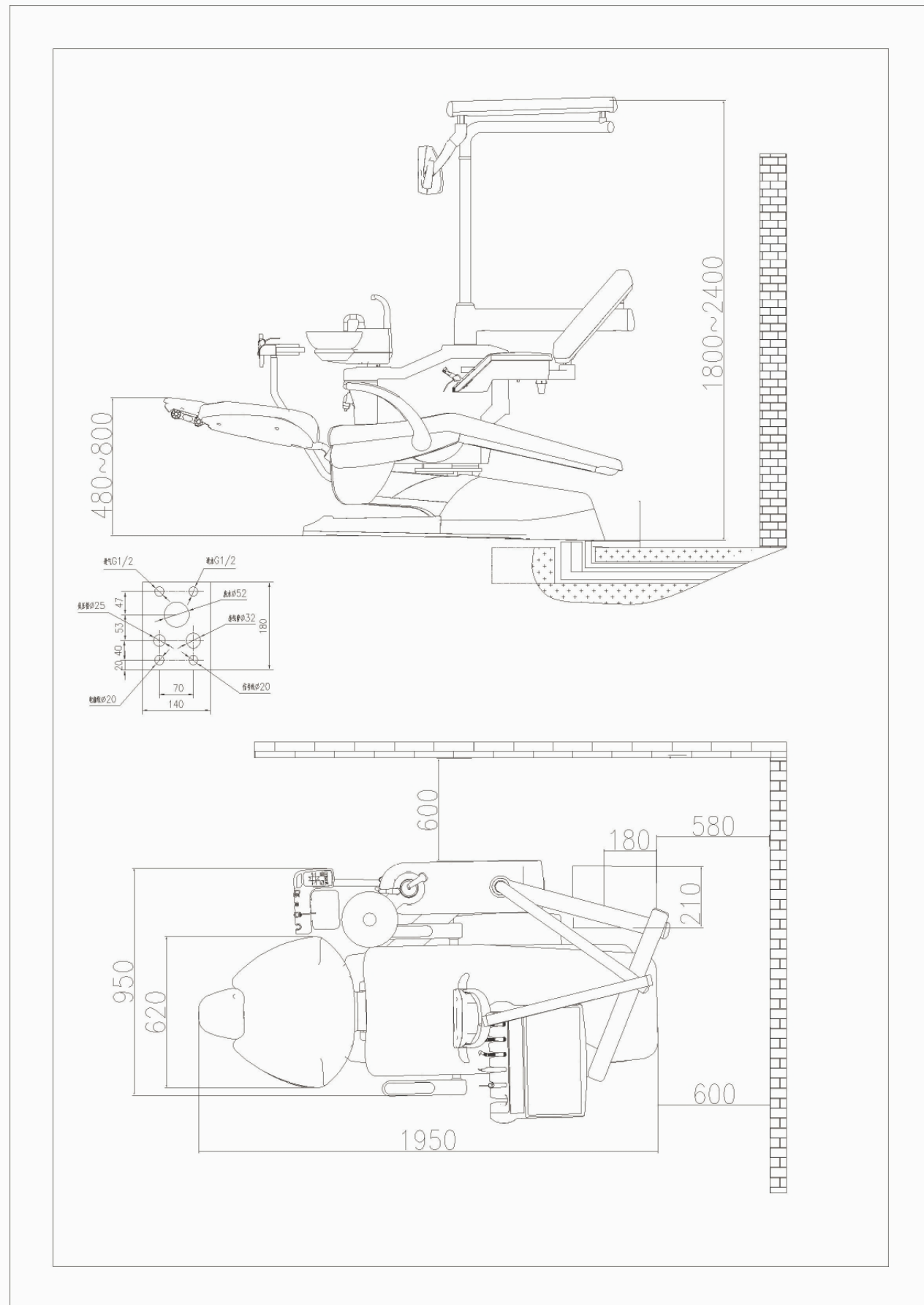

Standard configuration Micro fiber leather



PU for extra choices

